

GREEN&CLEAN INFECTION CONTROL

Hands

Surfaces

Special areas

Instruments



Safe

Clean

Efficient

Microorganisms are part of our daily lives. The risk of infection is usually not considered to be particularly high. In medical practices, however, the situation is different - the risk of infection is much higher here. The frequent change of patients and the nature of the work demand consistent hygiene measures. The responsibility of the physician towards his staff, the patient and himself calls for the use of highly effective cleaning and disinfecting agents.

With our GREEN&CLEAN products, METASYS offers you a comprehensive solution for all hygiene requirements and guidelines. Advanced active and conditioning agents ensure absolute safety and infection protection.

Ensure optimal hygiene in all areas of the dental practice with GREEN&CLEAN!

GREEN&CLEAN HD N GREEN&CLEAN HL | GREEN&CLEAN HC

GREEN&CLEAN HL - nourishing, antiseptic hand wash lotion

GREEN&CLEAN HD N - disinfectant for hygienic hand disinfection and surgical hand and forearm disinfection

GREEN&CLEAN HC - caring hand cream



GREEN&CLEAN WD


Alcoholic surface disinfection for small surfaces, such as those of straight & contra-angle handpieces


GREEN&CLEAN IK

Disinfection of spray water and spray air channels of turbines, straight & contra-angle handpieces

GREEN&CLEAN ID N

Disinfecting pre-cleaning of invasive and non-invasive instruments

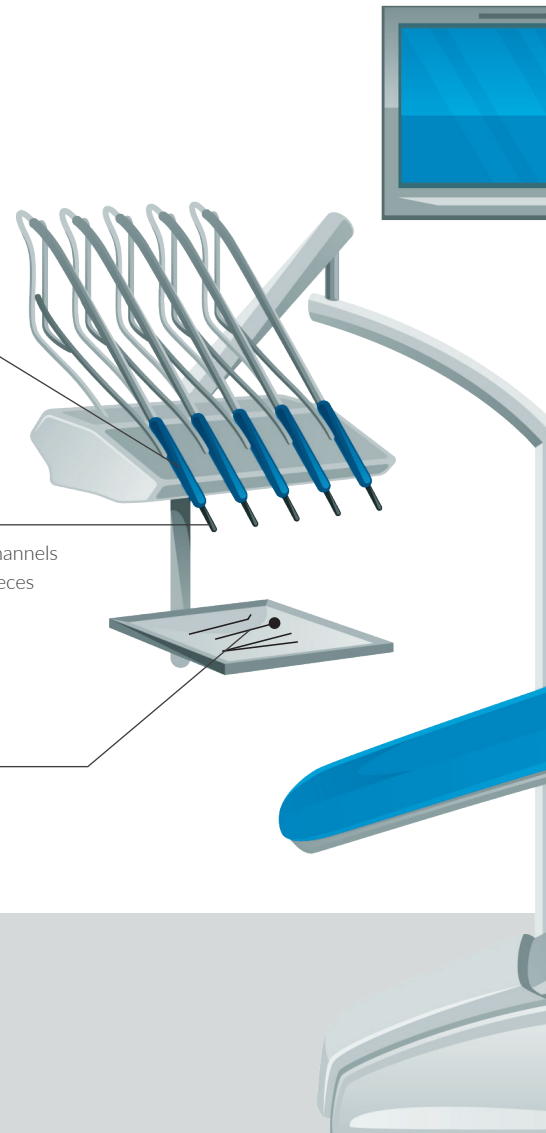

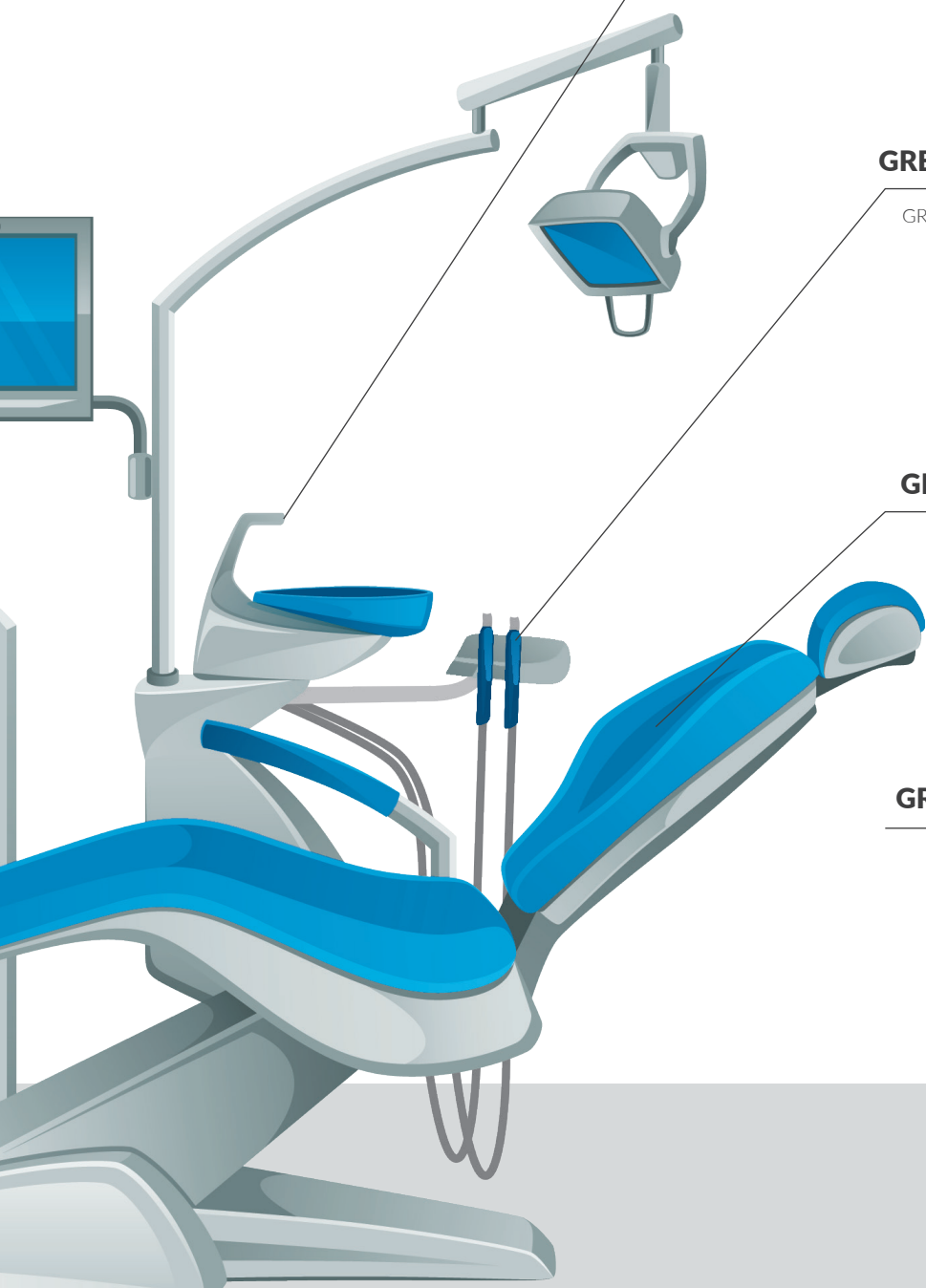

GREEN&CLEAN WK | GREEN&CLEAN BR

GREEN&CLEAN WK - decontamination and protection against lime deposits and biofilm in water-bearing lines

GREEN&CLEAN BR - biofilm removal

GREEN&CLEAN SD | GREEN&CLEAN MT

Alcoholic surface disinfection for alcohol-insensitive surfaces

GREEN&CLEAN M2 | GREEN&CLEAN CL N

GREEN&CLEAN M2 - cleaning, disinfection and deodorization of suction systems and amalgam separators

GREEN&CLEAN CL N - cleaning of suction systems

GREEN&CLEAN SK | GREEN&CLEAN MK

alcohol-free surface disinfection for alcohol-sensitive surfaces

GREEN&CLEAN MB | GREEN&CLEAN AD

GREEN&CLEAN MB - spittoon bowl cleaner

GREEN&CLEAN AD - impression disinfection

**GREEN
&CLEAN**

HD N

HAND DISINFECTION

Nourishing rub-in preparation for hygienic hand disinfection and surgical hand and forearm disinfection

- > high-quality alcohols prevent the skin from drying out
- > maintains the natural pH value of the skin
- > pleasant pine scent
- > VAH listed

Meets the following EN standards:

- > EN 1500 Hygienic hand disinfection
- > EN 12791 Surgical hand disinfection
- > EN 14476 Virucidal efficacy
- > EN 14348 Mycobactericidal efficacy

Modern active agent combination

Propan-1-ol, ethanol-1

Spectrum of efficacy

bactericidal, fungicidal (C. albicans), tuberculocidal, fully virucidal

Action times

Hygienic hand disinfection:

Rub 3 ml GREEN&CLEAN HD N on hands for 30 seconds.

Surgical hand disinfection:

Rub 10 ml GREEN&CLEAN HD N on hands and forearms for 1.5 minutes.

Order numbers

122000009 GREEN&CLEAN HD N, 6 x euro bottles 1000 ml

122000010 GREEN&CLEAN HD N, 6 x euro bottles 500 ml

*Use GREEN&CLEAN HD N with care. Always read label and product information before use.



DISPENSER SYSTEMS

Dispenser systems for hand disinfectants and liquid soaps



META TOUCHLESS

- > dispenser system for contactless, sensor-controlled dosing for 500 ml euro bottles
- > dosage: 1.5 ml per stroke
- > flexible positioning: As tabletop dispenser or wall-mounted
- > removable drip tray

Order number

121000019 META TOUCHLESS dispenser system for 500 ml bottles



META TOUCH

- > wall-mounted dispenser for 1000 ml euro bottles
- > dishwasher safe and autoclavable up to 134°C
- > dosage adjustable to 0.75 ml / 1.0 ml / 1.5 ml per stroke
- > housing made of sturdy aluminum
- > easy pump removal from the front

Order number

121000018 META TOUCH wall dispenser for 1000 ml bottles

**GREEN
&CLEAN**

HL

HYGIENIC HAND CLEANING

Nourishing hand washing lotion with broad antimicrobial effect

- > broad antiseptic effect
- > alkali-free, pH-neutral
- > protects the natural acid mantle of the skin
- > prevents the development of skin irritations
- > with allantoin as a skin care substance

Order numbers

122000013 1 x euro bottle 1000 ml

122000011 5 x bottles 750 ml + dosing bottle

122000012 6 x bottles 750 ml

122000014 1 x canister 5 l

121000001 Tap for canister 5 l



**GREEN
&CLEAN**

HC

HAND CARE

Caring hand cream with shea butter, ectoin and allantoin

- > suitable for sensitive skin
- > nourishes the skin and restores moisture, the skin looks supple and relaxed
- > counteracts skin cracks and brittleness
- > absorbs quickly

Order numbers

122000066 1 x 100 ml dispenser

122000007 6 x 100 ml dispenser



SURFACES

SURFACE DISINFECTION

After each treatment, surfaces that are near the patient and are contaminated by contact or aerosols, must be cleaned and disinfected. For aerosols, a dispersal radius of 2 meters must be taken into account. All surfaces farther away from the patient that have been visibly contaminated must also be disinfected immediately.

Once a working day, all other functional surfaces in the clinically used practice area must be treated with a surface disinfectant. For reasons of occupational safety, wipe disinfection is generally preferable to spray disinfection.



GREEN &CLEAN SK

SPRAY DISINFECTION FOR LARGE SURFACES

Alcohol-free spray disinfection for alcohol-sensitive surfaces of medical devices

- > alcohol-free
- > short action time
- > ready-to-use disinfection solution
- > especially suitable for alcohol-sensitive surfaces of medical devices
- > aldehyde-, phenol- and phosphate-free
- > CE 0044



Modern active agent combination

Quaternary ammonium compounds

| Spectrum of efficacy | Standard | Test information | Action time |
|--|----------------------|------------------|-------------|
| Bactericidal | EN 13727 | low load | 1 min |
| Levurocidal | EN 13624 | low load | 1 min |
| 4-field test (bactericidal, levurocidal) | EN 16615 | low load | 5 min |
| Limited virucidal | EN 14476 EN 16777 | low load | 1 min |

Order numbers GREEN&CLEAN SK

- 122000064 1 x spray bottle 750 ml, with foamer
- 122000050 6 x spray bottles 750 ml + 2 x foamer
- 122000051 1 x canister 5 l
- 121000001 Tap for canister 5 l

GREEN &CLEAN MK

WIPE DISINFECTION FOR LARGE SURFACES

Wipes for alcohol-free wipe disinfection of alcohol-sensitive surfaces of medical devices

- > extra large wipes (290 x 270 mm), for use with GREEN&CLEAN SK, for wipe disinfection of large surfaces of medical devices

Order numbers GREEN&CLEAN MK

- 122000044 1 x dispenser with 70 wipes + 5 x wipe rolls of 70 wipes each
1 x GREEN&CLEAN SK canister 5 l
- 122000042 2 x wipe rolls of 70 wipes each
- 122000043 1 x dispenser with 70 wipes + 5 x wipe rolls of 70 wipes each



SURFACES

GREEN & CLEAN SD

SPRAY DISINFECTION FOR LARGE SURFACES

Alcohol-based spray disinfection for alcohol-insensitive surfaces of medical devices

- > short action time
- > ready-to-use disinfection solution
- > formaldehyde-free
- > CE 0044

Modern active agent combination

Ethanol, 2-propanol, quaternary ammonium compounds, glyoxal



| Spectrum of efficacy | Standard | Test information | Action time |
|--|----------------------|------------------|-------------|
| Bactericidal | EN 13727 | high load | 15 s |
| Levurocidal | EN 13624 | high load | 15 s |
| Tuberculocidal | EN 14348 | high load | 1 min |
| 4-field test (bactericidal, levurocidal) | EN 16615 | high load | 1 min |
| Limited virucidal | EN 14476 EN 16777 | low load | 1 min |

Order numbers GREEN&CLEAN SD

- 122000047 5 x bottles 750 ml + 2 x spray heads
- 122000045 6 x bottles 750 ml
- 122000048 1 x canister 5 l
- 122000046 2 x bottles 750 ml + 1 x spray head
- 121000001 Tap for canister 5 ml

GREEN & CLEAN MT

WIPE DISINFECTION FOR LARGE SURFACES

Alcohol-based wipe disinfection for alcohol-insensitive surfaces of medical devices

- > extra large wipes (290 x 270 mm), for use with GREEN&CLEAN SD, for wipe disinfection of large surfaces of medical devices

Order numbers GREEN&CLEAN MT

- 122000044 1 x dispenser with 70 wipes + 5 x wipe rolls with 70 wipes each
- 1 x 5 l GREEN&CLEAN SD canister
- 122000042 2 x wipe rolls of 70 wipes each
- 122000043 1 x dispenser with 70 wipes + 5 x wipe rolls of 70 wipes each



GREEN
& CLEAN

WD

WIPE DISINFECTION FOR SMALL SURFACES

Alcohol-based disinfection wipes made of tear-resistant fleece for disinfection of surfaces of medical devices (e.g. outer surfaces of handpieces & contra-angles)

- > short action time
 - > formaldehyde-free
 - > with hygienic safety seal
- > CE 0044

Modern active agent combination

Ethanol, 2-propanol, quaternary ammonium compounds, glyoxal

| Spectrum of efficacy | Standard | Test information | Action time |
|--|----------------------|------------------|-------------|
| Bactericidal | EN 13727 | high load | 30 s |
| Levurocidal | EN 13624 | high load | 30 s |
| Tuberculocidal | EN 14348 | high load | 30 s |
| 4-field test (bactericidal, levurocidal) | EN 16615 | high load | 1 min |
| Limited virucidal | EN 14476 EN 16777 | high load | 30 s |

Order numbers

- 122000052 1 x dispenser with 120 wipes
- 122000053 1 x dispenser with 120 wipes + 3 x bags of 120 wipes each
- 122000055 6 x bags of 120 wipes each
- 121000012 GREEN&CLEAN WD dispenser, empty



**GREEN
& CLEAN**

M2

SUCTION SYSTEM DISINFECTION

Dual, active defoaming concentrate for the daily cleaning, disinfection and deodorization of suction systems and amalgam separators

- > dual system (red / green)
- > excellent active defoaming
- > easy handling through well-designed dispenser
- > enzymatic product with excellent protein dissolving power
- > ideal for all amalgam separators and water-air separations
- > CE 0044



Modern active agent combination

Alkylamine, quaternary ammonium compounds, enzymes

| GREEN&CLEAN M2 green Spectrum of efficacy | Standard | Test information | Action time | solution |
|---|----------|------------------|-------------|-------------|
| Bactericidal | EN 13727 | high load | 10 min | 1% solution |
| Levurocidal | EN 13624 | high load | 10 min | 1% solution |
| Tuberculocidal | EN 14348 | high load | 240 min | 1% solution |
| Germ carrier test - bactericidal | EN 14561 | high load | 15 min | 1% solution |
| Germ carrier test - levurocidal | EN 14562 | high load | 15 min | 1% solution |
| Limited virucidal | EN 14476 | high load | 60 min | 1% solution |

| GREEN&CLEAN M2 red Spectrum of efficacy | Standard | Test information | Action time | solution |
|---|----------|------------------|-------------|-------------|
| Bactericidal | EN 13727 | high load | 10 min | 1% solution |
| Levurocidal | EN 13624 | high load | 10 min | 1% solution |
| Tuberculocidal | EN 14348 | high load | 240 min | 1% solution |
| Germ carrier test - bactericidal | EN 14561 | high load | 15 min | 1% solution |
| Germ carrier test - levurocidal | EN 14562 | high load | 15 min | 1% solution |
| Limited virucidal | EN 14476 | high load | 60 min | 1% solution |

Order numbers

| | |
|-----------|---|
| 122000026 | 2 x bottles 500 ml (1 x red, 1 x green) + dispenser |
| 122000027 | 4 x bottles 500 ml (2 x red, 2 x green) |
| 122000028 | 2 x bottles 500 ml (1 x red, 1 x green) |
| 122000031 | GREEN&CLEAN M2 red, 25 x bottles 500 ml |
| 122000030 | GREEN&CLEAN M2 green, 25 x bottles 500 ml |
| 121000009 | Dispenser for GREEN&CLEAN M2 |

SIMPLE - EFFICIENT - UNIQUE

Disinfection and cleaning of the suction system with GREEN&CLEAN M2

- > Efficient: Broad spectrum of efficacy and high cleaning power ensure profound hygiene
- > Simple: Ergonomic dosing dispenser enables rapid disinfection process
- > Unique: Distinctive color combination as an expression of different, active ingredients for long-lasting safety



SPECIAL AREAS

GREEN
& CLEAN

CL N

SUCTION SYSTEM CLEANING

Concentrate for the daily cleaning of suction systems

- > improved defoaming due to modern combination of active ingredients
- > excellent protein dissolving power due to special enzyme mixture
- > pleasant scent in the treatment rooms
- > simple dosing thanks to double-neck bottle
- > pH-neutral user solution
- > CE

Modern active agent combination

Alkyl ether carboxylic acid, subtilisin

Order numbers

122000004 2 x bottles 1000 ml

121000005 Mixing container for GREEN&CLEAN CL N



GREEN & CLEAN AD

IMPRESSION DISINFECTION

Ready-to-use solution for the sustained immersion disinfection of all prosthetics, worn dentures and dental impressions

- > short action time
- > suitable for alginates, silicones, polyethers and polysulfide
- > high material compatibility
- > formaldehyde and phenol free
- > standing time: 7 days
- > CE 0044

Modern active agent combination

Quaternary ammonium compounds, alkylamine, tensides



| Active agent | Test germ | Standard | Test information | Action time |
|----------------------------------|------------------------------------|----------|------------------------|-------------|
| Bactericidal | S. aureus, E. hirae, P. aeruginosa | EN 13727 | high load ¹ | 2 min |
| Levurocidal | C. albicans | EN 13624 | high load ¹ | 2 min |
| Germ carrier test - bactericidal | S. aureus, E. hirae, P. aeruginosa | EN 14561 | high load ¹ | 5 min |
| Germ carrier test - levurocidal | C. albicans | EN 14562 | high load ¹ | 5 min |
| Limited virucidal | Vacciniavirus | EN 14476 | low load ² | 15 sec |
| | | EN 16777 | low load ² | 1 min |

¹ Bovine albumin 3g/l + sheep erythrocytes 3ml/l

² Bovine albumin 0.3g/l

Order numbers

- 122000001 1 x canister 5 l
- 121000003 Disinfection tray 1.25 l with cover and strainer
- 121000001 Tap for 5 l canister

GREEN & CLEAN MB

SPITTOON BOWL CLEANING

Spittoon bowl cleaner based on citric acid and enzymes, with integrated brush

- > excellent protein dissolving power due to enzymatic ingredients
- > with special defoamer
- > CE

Modern active agent combination

Quaternary ammonium compounds, tensides, enzymes

Order numbers

- 122000035 Dispenser bottle 180 ml, with brush
- 122000032 2 x bottles 500 ml + 2 empty dispenser bottles + 2 brushes
- 122000033 4 x bottles 500 ml + 4 brushes
- 122000034 4 x bottles 500 ml
- 122000037 4 x spare brushes
- 122000036 Dispenser bottle with brush, empty



SPECIAL AREAS

GREEN
& CLEAN

BR

BIOFILM REMOVAL

4% solution for the removal of biofilm from the water-bearing lines of the dental unit with the BR biofilm cleaning device or bottle systems

- > effective biofilm removal in only 30 minutes
- > recommended for the refurbishment of dental units as well as before the installation of a water disinfection system
- > simple application due to colored solution
- > CE 0044

Active agent

Hydrogen peroxide

Order number

122000002 Biofilm cleaning solution, 2 x 1000 ml bottles



GREEN
& CLEAN

H1

CLEANING AGENT FOR THE H1 HYGIENE SYSTEM

Decontamination concentrate for the automatic METASYS hygiene system H1

- > for the continuous cleaning of the entire suction system
- > active defoaming
- > CE

Modern active agent combination

Quaternary ammonium compounds, special tensides, enzymes

Order number

122000005 4 x bags 500 ml



METASAN®

P

WATER DISINFECTION

1.41% solution based on hydrogen peroxide, suitable for disinfection of the treatment water in Dentsply-Sirona treatment units

- > continuous disinfection of the treatment water and the water-carrying pipes
- > continuous decontamination of the service water and the water-bearing lines
- > broad spectrum of efficacy
- > CE 0044

Active agent

Hydrogen peroxide

Order numbers

122000059 METASAN® P 6 x bottles 1000 ml



**GREEN
&CLEAN**

WK

WATER DECONTAMINATION

2% solution for decontamination and protection against lime deposits in the water-bearing lines in a dental unit based on hydrogen peroxide

- > wide spectrum of action
- > prevention of lime deposits
- > also suitable for bottle systems
- > CE 0044

Active agent

Hydrogen peroxide

GREEN&CLEAN WK test strips

- > Test strips for the visual inspection of the functionality and effectiveness of the water decontamination device
- > the color scale on the test strip indicates the concentration of GREEN&CLEAN WK in the process water

dosing aid for bottle systems

- > suitable for GREEN&CLEAN WK 750 ml

Order numbers

122000057 4 x bottles 750 ml
122000056 6 x bottles 1000 ml
121000015 WK test strips, 25 pcs.

121000017 WK test strips, 50 pcs.
121000014 Dosing aid for bottle systems
(for 750 ml bottle)



**GREEN
& CLEAN**

IK

DISINFECTION OF STRAIGHT & CONTRA-ANGLE HANDPIECES

Ready-to-use, alcohol-based disinfectant spray for the spray water and spray air channels of turbines, straight & contra-angle handpieces

- > short action time
- > increased application safety due to fixed adapters
- > one spray can is sufficient for approx. 300 applications at approx. 2 seconds spraying time
- > one universal adapter for all straight and contra-angle couplings and own turbine couplings
- > does not replace oil spray recommended by the manufacturer!
- > CE 0044



Modern active agent combination

1-propanol, 2-propanol, 2-butanol, quaternary ammonium compounds

| Spectrum of efficacy | Standard | Test information | Action time |
|------------------------------------|----------------------|------------------|-------------|
| Bactericidal | EN 13727 | high load | 1 min |
| Levurocidal | EN 13624 | high load | 1 min |
| Levurocidal | EN 14348 | high load | 1 min |
| Germ carrier test - bactericidal | EN 14561 | high load | 2 min |
| Germ carrier test - levurocidal | EN 14562 | high load | 1 min |
| Germ carrier test - tuberculocidal | EN14563 | high load | 1 min |
| Limited virucidal | EN 14476 EN 16777 | high load | 30 s |

Order numbers

2 x 200 ml spray cans incl. pre-assembled adapter:

- | | |
|-----------|--|
| 122000017 | Hand pieces & contra-angles (universal) |
| 122000025 | Turbine coupling W&H |
| 122000021 | Turbine coupling Bien Air |
| 122000023 | Turbine coupling KaVo |
| 122000024 | Turbine coupling Dentsply Sirona |
| 122000022 | Turbine coupling Castellini |
| 122000018 | Turbine coupling NSK |
| 122000020 | Handpieces & contra-angles (universal) Dentsply Sirona T1 Classic |
| 122000019 | Turbine coupling NSK Quick |



GREEN
&CLEAN

IDN

INSTRUMENT DISINFECTION

Highly effective concentrate for the disinfecting pre-cleaning of invasive and non-invasive instruments

- > 5 l concentrate yields up to 500 l working solution at 1% concentration
- > aldehyde and phenol free
- > VAH listed
- > standing time: 7 days
- > CE 0044

Modern active agent combination

Quaternary ammonium compounds, alkylamine, corrosion inhibitors



| Spectrum of efficacy | Standard | Test information | Action time | solution |
|------------------------------------|----------|------------------|-------------|---------------|
| Bactericidal | EN 13727 | high load | 5 min | 2% solution |
| Levurocidal | EN 13624 | high load | 5 min | 1% solution |
| Tuberculocidal | EN 14348 | high load | 60 min | 1% solution |
| Germ carrier test - bactericidal | EN 14561 | high load | 15 min | 1% solution |
| Germ carrier test - levurocidal | EN 14562 | high load | 15 min | 1% solution |
| Germ carrier test - tuberculocidal | EN 14563 | high load | 60 min | 0.5% solution |
| Limited virucidal | EN 14476 | low load | 15 min | 0.5% solution |

Order numbers

| | | | |
|-----------|---------------------------------------|-----------|--------------------------------------|
| 122000015 | 4 x bottles 500 ml + measuring beaker | 121000007 | Measuring beaker for GREEN&CLEAN IDN |
| 122000016 | 1 x canister 5 l | 121000001 | Tap for canister 5 l |



GREEN&CLEAN



HANDS

| | | Cleaning | Disinfection | Alcoholic | Ready for use |
|------------------|------------------------|----------|--------------|-----------|---------------|
| GREEN&CLEAN HL | Hygienic hand cleaning | ■ | | | ■ |
| GREEN&CLEAN HD N | Hand disinfection | | ■ | ■ | ■ |
| GREEN&CLEAN HC | Hand care | | | | ■ |

SURFACES

| | | | | | |
|----------------|--------------------|---|---|---|---|
| GREEN&CLEAN SK | Spray disinfection | ■ | ■ | | ■ |
| GREEN&CLEAN MK | Wipe disinfection | ■ | ■ | | |
| GREEN&CLEAN WD | Wipe disinfection | ■ | ■ | ■ | ■ |
| GREEN&CLEAN SD | Spray disinfection | ■ | ■ | ■ | ■ |
| GREEN&CLEAN MT | Wipe disinfection | ■ | ■ | ■ | |

SPECIAL AREAS

| | | | | | |
|----------------------|---|---|---|--|---|
| GREEN&CLEAN M2 green | Suction system disinfection | ■ | ■ | | |
| GREEN&CLEAN M2 red | Suction system disinfection | ■ | ■ | | |
| GREEN&CLEAN H1 | Cleaning agent for H1 | ■ | | | ■ |
| GREEN&CLEAN CL N | Suction system cleaning | ■ | | | |
| GREEN&CLEAN MB | Spittoon bowl cleaning | ■ | | | ■ |
| GREEN&CLEAN WK | Water disinfection and biofilm prevention | | ■ | | |
| METASAN® P | Water disinfection and biofilm prevention | | ■ | | ■ |
| GREEN&CLEAN BR | Biofilm removal | ■ | ■ | | ■ |
| GREEN&CLEAN AD | Impression disinfection | ■ | ■ | | ■ |

INSTRUMENTS

| | | | | | |
|------------------|--|---|---|---|---|
| GREEN&CLEAN IK | Disinfection of straight & contra-angle handpieces, turbines | ■ | ■ | ■ | ■ |
| GREEN&CLEAN ID N | Instrument disinfection | ■ | ■ | | |



Tested according to the following EN standards:

| EN 13623 | EN 13624 | EN 13727 | EN 14348 | EN 14476 | EN 14561 | EN 14562 | EN 14563 | EN 16615 | EN 16777 | Available as |
|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|--------------------------|
| | | | | | | | | | | 750 ml, 1000 ml, 5000 ml |
| | | | | | | | | | | 500 ml, 1000 ml |
| | | | | | | | | | | 100 ml |
| | ■ | ■ | | ■ | | | | ■ | ■ | 750 ml, 5000 ml |
| | | | | | | | | | | 70 wipes |
| | ■ | ■ | ■ | ■ | | | | ■ | ■ | 120 wipes |
| | ■ | ■ | ■ | ■ | | | | ■ | ■ | 750 ml, 5000 ml |
| | | | | | | | | | | 70 wipes |
| | ■ | ■ | ■ | ■ | ■ | ■ | | | | 500 ml |
| | ■ | ■ | ■ | ■ | ■ | ■ | | | | 500 ml |
| | | | | | | | | | | 500 ml |
| | | | | | | | | | | 1000 ml |
| | | | | | | | | | | 180 ml, 500 ml |
| ■ | | | | | | | | | | 750 ml, 1000 ml |
| ■ | | | | | | | | | | 1000 ml |
| ■ | | | | | | | | | | 1000 ml |
| | ■ | ■ | | ■ | ■ | ■ | | | ■ | 5000 ml |
| | ■ | ■ | ■ | ■ | ■ | ■ | ■ | | ■ | 200 ml |
| | ■ | ■ | ■ | ■ | ■ | ■ | ■ | | | 500 ml, 5000 ml |

APPLICATION AIDS



Disinfection tray 1.25 l

Order no. 121000003



Bur disinfection box 200 ml

Order no. 121000002



META TOUCH

Order no. 121000018



META TOUCHLESS

Order no. 121000019


Tap for canister 5 l

Order no. 121000001

GREEN&CLEAN WK test strips

Order no. 121000015 - WK test strips, half-year requirement, 25 pcs.

Order no. 121000017 - WK test strips, annual requirement, 50 pcs.

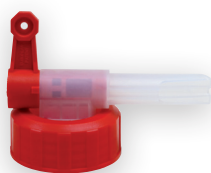

Dosing aid for bottle systems

Order no. 121000014



 **METASYS**
protect what you need

METASYS Medizintechnik GmbH

METASYS logistics & collection GmbH

METASYS recycling & processing GmbH

Florianistraße 3, 6063 Rum bei Innsbruck, Austria
T +43 512 205420

www.metasys.com | info@metasys.com

